

PRACTICAL

RUGGED

SMART

LABQUEST  
BY BOROSIL®

## DATA LOGGER PRINT PORTAL



Manufactured in  
**BOROSIL**  
TECHNOLOGIES



DATA LOGGER



INTELLIGENT



PRECISION






# Borosil Data Logger

Data logging is the recording of collected data over a period of time.

Borosil instruments with USB connectivity e.g. Industrial Stirrers and premium Water Baths which facilitate data logging.

Borosil Data Logger continuously monitors and records the activity performed by connected instruments over a period of time. Users can take the print out for specified period and operation. PNP010 a reliable, cost and time-efficient monitoring solution for any measuring & recording progress of set parameters offering key advantages by providing high accuracy and ease of use.

- Display 128x64 GLCD, 2 USB ports, internal memory 16 GB, alphanumeric keypad, RTC with battery life up to 5 yrs
- Login feature supports administrator and up to 99 users
- Data logger capacity up to 32 million files (@file size...4KB)
- Print paper size width 48mm with maximum 32 characters
- Power supply 12V DC adapter as a power source
- Any time print option for logged data records
- Customisable print format with company name, lab/department name and process parameters
- 12V DC adapter, battery consist of 4400 mAh battery with battery backup of continuous 5 hours
- Users can assign instrument ID (Connected Instrument ID) from PNP010 device

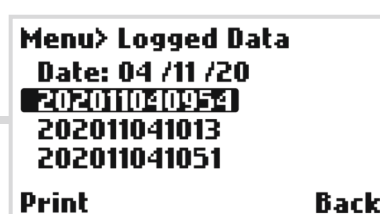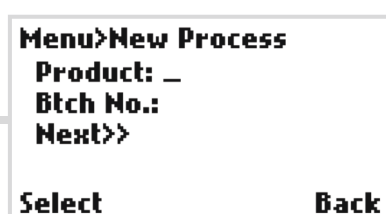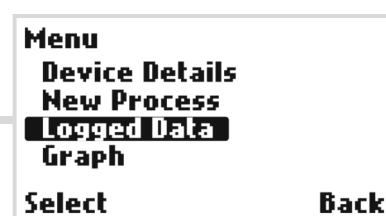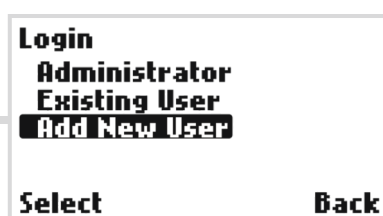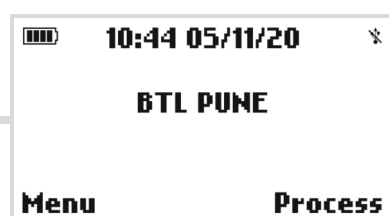

**BOROSIL®**  
BOROSIL LIMITED

1101, Crescenzo, G-Block, Opp. MCA Club,  
Bandra Kurla Complex, Bandra (E), Mumbai - 400 051, India  
T : +91 22 6740 6300 F : +91 22 6740 6514  
E : borosil@borosil.com W : www.borosil.com

Manufactured by :  
**BOROSIL TECHNOLOGIES LTD.**  
21-24, Plot No. 8, Indialand Global Industrial Park,  
Hinjewadi Phase - 1, Pune - 411 057