



Hybrid e-Dispenser

The Ultimate Laboratory Dispensing solution

Top-Notch Features for Precision and Accuracy

Precision Control Knob for effortlessly toggling between 3 preset speeds.


Distinct Modes for Dispensing and Titration.

Accuracy and Precision Complying with ISO 8655

USB interface facilitates transfer of data to PC

QR scanner for easy access to Chemical compatibility chart

Includes a range of adapters for mounting on diverse bottle neck

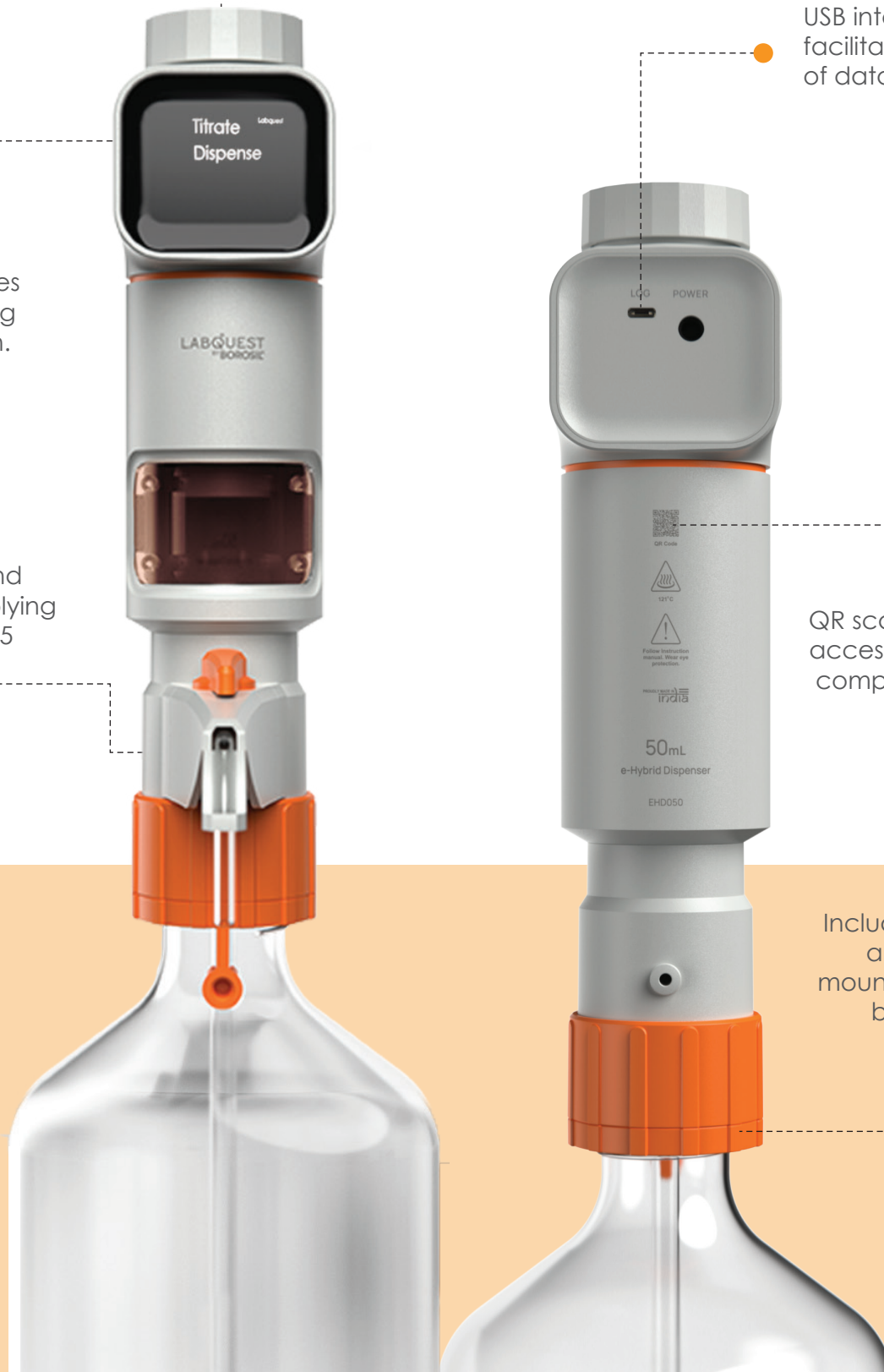

SMART

RUGGED

PRACTICAL



Select between the 3 preset speed- rapid dispensing , intermediate dispensing and meticulous dispensing for precise control



Amber colored and transparent Inspection window



Enhanced Visibility - Large and clear display ensuring visible reading from every angle.

Reliable results with no downtime with inbuilt calibration feature

SMART

RUGGED

PRACTICAL

Technical Specifications

Volume	50mL
Volume Accuracy in Dispense mode	+/- 0.30ml
Volume Accuracy in Titration mode	+/- 0.05ml
Dispense Incremental Value	0.2mL
Dispenser Modes	Time based, Fixed volume dispensing
Speed for Titration	3 Calibrated Pre-set Speeds
Operation Mode	Dispensing + Titration
Viewing Window	Amber & Clear Transparent
Operating method	Piston Operated
Inbuilt Calibration	Yes
Adapters	GL28, GL32, GL36, GL38, GL40 & GL45
Connectivity for data Storage	USB interface for PC connectivity, enables storage of upto 10 readings



BOROSIL®
Scientific

1101, Crescenzo, G-Block, Opp. MCA Club,
Bandra Kurla Complex, Bandra (E), Mumbai - 400 051, India.
T +91 22 6740 6300 F +91 22 6740 6514
E bsl@borosil.com W www.borosilscientific.com

PROUDLY MADE IN 
india



SMART

RUGGED

PRACTICAL